IMPORTANT SAFETY INSTRUCTIONS

When using electrical appliances, especially when children are present, basic safety precautions should always be taken, including the following:

READ ALL INSTRUCTIONS BEFORE USING SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL

DANGER — Any appliance is electrically live, even when the switch is off.

WARNING — To reduce the risk of burns, electrocution, fire, injury or death to persons:

1. Always unplug appliance from electrical outlet immediately after using (except when charging tool). Do not allow to be used as a toy.
2. Do not place or store appliance where it could fall or be pulled into a tub or sink.
3. Do not use appliance while bathing or in the shower.
4. Do not place appliance in or drop into water or other liquids.
5. If appliance has fallen into water, do not reach for it. Unplug immediately.
6. Unplug this appliance before cleaning.
7. Do not use device if it is not working properly or if it has been dropped, damaged, or dropped into water.
8. Do not place or store the charging cable where it can fall or be pulled into a tub or sink.
9. Do not immerse charging cable in water or any other liquids.
10. If charging cable falls into water or another liquid,

11. Not for use by children.
12. Do not pull, twist or wrap the cord around the appliance.
13. Never drop or insert any object into any opening on the appliance.
14. THIS PRODUCT IS FOR HOUSEHOLD USE ONLY.
15. Do not leave appliance unattended when plugged in. Unplug from outlet when not in use and before cleaning and servicing.
16. Do not unplug by pulling on cord. To unplug, grasp the plug, not the cord.
17. Do not handle plug or appliance with wet hands.

CAUTION —

1. Dry location use only. Do not expose to liquids, vapor, or rain.
2. Ensure jack is plugged into the tool before plugging the USB cable into the outlet.
3. Do not operate appliance with a damaged housing or use charging cord that has been damaged.
4. Always unplug the unit before cleaning or when not in use.
5. Be sure voltage is correct before plugging in. The rated input is 5V DC 1.0A.

The SKN by conair pore purifier advanced™ micro-dermabrasion tool handle, attachment heads and charging cable have no serviceable parts. If any component is damaged, it must be replaced. Call 1-800-3-CONAIR to purchase replacements.

SAVE THESE INSTRUCTIONS

SENSITIVITY AND CAUTIONS SPECIAL CAUTIONS FOR DELICATE AREAS

The SKN by conair™ Microdermabrasion Tool is designed to gently clean your pores and exfoliate your skin without irritation, for healthy-looking results.


If you are under a doctor's care for any skin condition, or if you think you have very sensitive skin, please consult your physician before use.

We recommend starting on the LOW SPEED/SUCTION setting and increasing to HIGH SPEED/SUCTION, depending on how your skin reacts to treatment. If rash or irritation occurs, stop use.

This is not a medical device and is not intended to be used to treat any medical conditions.

We recommend not sharing the device. If you experience any pain or discomfort during use, stop using and consult your physician.

SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL

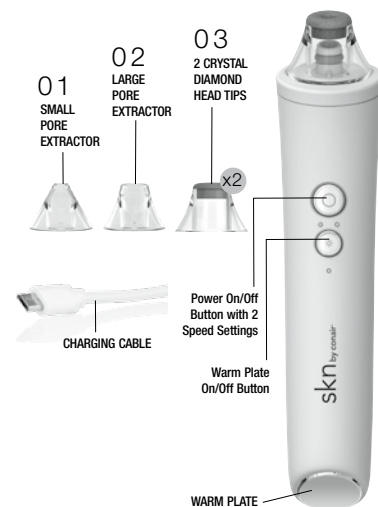USER MAINTENANCE INSTRUCTIONS

SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL

It is recommended to clean attachment heads and device after each use. If not clean, performance could be affected. Change the filter if the suction strength is weak. Remove attachment head, peel off metal filter, stick on new filter, and reattach the microdermabrasion tip.

1. After each use, wipe the handle with a damp cloth. The attachment heads and filters can be rinsed with soapy water, but the handle is not waterproof and cannot be rinsed under water or any other liquid. Also, do not put unit in dishwasher for cleaning.

2. After cleaning, wipe all parts with a dry cloth and store in a clean, dry place.

TROUBLE-SHOOTING GUIDE

SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL IS NOT CHARGING NO WHITE FLASHING CHARGING LIGHT.	<ol style="list-style-type: none"> 1. Make sure the contacts on the unit and charging cable are clean. See User Maintenance Instructions for cleaning procedure. 2. Make sure the unit is plugged in properly. 3. Check the outlet and make sure the power is on. Most bathroom outlets have a test and reset button. Press the reset button and recheck the outlet.
THE CHARGING LIGHT DOES NOT CHANGE TO SOLID WHITE.	
SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL WILL NOT TURN OFF.	You must hold the top power button for more than 1 second to turn the unit off.
SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL OR ATTACHMENTS FALL OFF DURING USE.	Make sure the attachments are snapped on properly.

BATTERY DISPOSAL

The battery is sealed inside the SKN by conair pore purifier advanced™ microdermabrasion tool so it cannot be easily removed for recycling. The entire SKN by conair pore purifier advanced™ microdermabrasion tool should be returned when recycling the battery for this product.

Bring the device to the nearest collection center for proper disposal. DO NOT throw the device into the general municipal waste system.

CONTAINS LITHIUM-ION BATTERY. MUST BE COLLECTED, RECYCLED OR DISPOSED OF PROPERLY. DO NOT DISPOSE OF IN THE MUNICIPAL WASTE.

RBRC™ is a trademark of the Rechargeable Battery Recycling Corporation.

To find a battery recycling center near you, call 1-800-822-8837 or go to www.call2recycle.org/locator/.

FEATURES

TWO SPEED OPTIONS

The two speeds available on your SKN by conair pore purifier advanced™ microdermabrasion tool allow you to customize your skincare routine to suit your individual needs.

WARM PLATE

The warm plate on the bottom of the tool helps to open pores, leaving skin feeling bright and smooth.

ATTACHMENTS

SMALL PORE EXTRACTOR: To remove blackheads & whiteheads.

LARGE PORE EXTRACTOR: To target clogged pores.
2 CRYSTAL DIAMOND HEAD TIPS: For deep exfoliation.

NOTE: Your SKN by Conair may only include certain attachments. See retailer or call 1-800-3-CONAIR for additional attachments.

RECHARGEABLE

The SKN by conair pore purifier advanced™ microdermabrasion tool is cordless and requires a full charge before use. After tool is fully charged, it can provide up to 120 minutes of use without recharging.

OPERATING PRECAUTIONS

This device has a built-in rechargeable lithium-ion battery. Do not throw device into fire, apply heat to it, or charge, use in a high-temperature environment.

CHARGING THE SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL

1. Before using the SKN by conair pore purifier advanced™ microdermabrasion tool for the first time, it must be fully charged for 2.5–3.5 hours. Plug the DC jack into the side of the unit, then plug power cable and adapter into wall outlet. The white LED will begin to flash, indicating that the tool is charging.

2. The charging indicator will FLASH WHITE when charging. Once fully charged, the light will illuminate SOLID WHITE. Always check to be sure that the charging light illuminates when you plug the charging cable into the outlet.

3. The white light will start flashing for 5 seconds and then shut off when the unit reaches low battery capacity during operation.

4. The adapter may become warm while charging — this is normal. Once fully charged, the device can provide 120 minutes of use.

5. In order to maintain the optimum capacity of the rechargeable battery, the SKN by conair pore purifier advanced™ microdermabrasion tool should be given a full 2.5–3.5-hour charge every 3 months.

6. Remove charging cable from unit and unplug adapter from outlet when not in use. Do not unplug the adapter from the wall when it is still connected to the unit. This will cause the battery to discharge.

STARTING THE SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL AND SELECTING THE SPEED/SUCTION LEVEL

Press the power button once to turn on the device in low speed/suction. For a deeper exfoliation or stronger suction, press button 2 times for high speed/suction. The corresponding light will illuminate.

STOPPING THE SKN BY CONAIR PORE PURIFIER ADVANCED™ MICRO-DERMABRASION TOOL

At any time, the SKN by conair pore purifier advanced™ microdermabrasion tool can be stopped by pressing and holding the power switch button for more than 1 second.

TURNING ON THE WARM PLATE

Press and hold the Warm Plate Button for 1 second to activate the WARM PLATE. After 5 minutes in this setting, the device will go into standby mode.

TURNING OFF THE WARM PLATE

At any time, the WARM PLATE can be stopped by pressing and holding the Warm Plate Button for more than 1 second.

ATTACHING AND REMOVING MICRO-DERMABRASION/PORE EXTRACTOR HEADS

To attach the different attachments, hold the unit and push the attachment onto the post until it snaps in place. To remove attachment head, simply pull attachment off of the unit.

FCC MARKINGS

WARNING: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

CAUTION: Any changes or modifications not expressly approved by Conair LLC, the manufacturer of this product and the party responsible for compliance, could void the user's authority to operate the equipment.

OPERATING INSTRUCTIONS

Warm Plate and Attachments

1. Thoroughly cleanse and dry the skin.
2. To open pores for better results, use the warm plate before using the microdermabrasion or pore extractor attachments.
3. Test warm plate on a small area of the neck or forearm.
4. Applying light pressure, slide the warm plate gently over the facial skin. Apply in such a way that it gently warms the skin.
5. When warm treatment is complete, turn off and continue with microdermabrasion attachments.
6. Test microdermabrasion attachment on a small area of the neck or forearm.
7. Pull the skin taut with your free hand to allow the device to smoothly glide over the skin.
8. Move device in an upward, fluid motion across the skin. Do not hover in any spots.
9. Immediately after treatment wash and dry the skin.
10. Apply a calming toner and recovery moisturizer.
11. Clean the attachment heads and unit, including the warm plate, after each use. See User Maintenance Instructions for more information on cleaning the device and heads.
12. Wait 6 to 7 days between treatments.

NOTE: Do not do more than two passes over the same area during treatment.

SUPPLIER'S DECLARATION OF CONFORMITY

Brand: SKN by Conair
Model No.: MD03
Description: Pore Purifier Advanced™ Microdermabrasion Tool
Responsible Party: Conair LLC, 1 Cummings Point Road, Stamford, CT 06902 (203) 351-9000
Standards: FCC Part 15 Subpart B, ANSI C63.4-2014

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesirable operation.

DATE OF ISSUE: T8D

Please register this product at register.conair.com

SERVICE CENTER CONAIR LLC

Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

Visit us on the web at www.conair.com

©2022 Conair LLC
Glendale, AZ 85307
stock.image

MD03 Models
22PS082977

IB-17605

skin
by conair™

DEEP CLEAN +
EXFOLIATE

pore purifier advanced™
microdermabrasion tool

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Parte responsable: Conair LLC, 1 Cummings Point Road, Stamford, CT 06902 (203) 351-9000